

Product Parameter

Front wheel	8" wheel
Rear Wheel	12.5" wheel
Net Weight	26.5kg
Brake	Electromagnetic Braking system
Batteries	Lithium Battery 24V 11Ah
Motors (pair)	Brushless Motor 180W
Charger	24V 2A/3A/5A
Controller	Reliable Brushless Controller
Maximum Range	12-15 km
Reclining	Yes (4 positions adjustment)
Ground Clearance	6 cm
Max Speed	7 km/h
Charging Hour	5-6 Hours
Max Slope	12 degrees
Weight Capacity	125 kg
Dimension (Folded)	79*37*74cm
Dimension (Unfold)	109*63.5*98cm

1. Range and speed vary with user weight terrain type, battery charge battery condition and tyre pressure.
2. Due to manufacturing tolerances and continual product improvement, this specification can be subject to variance of + or - 3%.
3. Battery weight may vary based on manufacturer.

TO BE A YATTL AGENT

2020

Requirements

1. At least 1 lawful and registered local shop or sale place managed by yourself.
2. Legally registered company.
3. Have some knowledge and strong interest in electric wheelchairs.

Our Support

1. Agency area protection.
2. The best quality and competitive price.
3. Prior rights to the new products with new opportunity.
4. No minimum ordering quantity.
5. Soonest response on any of your questions.
6. Free parts supply support policy.

JOIN US



YATTL INDUSTRY CO., LTD.



Tel: +86 755 8170 7199

Web: www.yattll.com

Email: info@yattll.com

Factory: No.1 Lingchang Road, PingShan Village, TangXia Town, DongGuang City, GuangDong, China.

Office: RM.1403, BLK.B of Duocai Office Complex, Gu-anlan High-tech Zone, Longhua DIST., Shenzhen, China.

New product



Latest reclining power wheelchair

Model: **YEC35A**



www.yattll.com |    

Detail display chart



NEW PRODUCT

YEC35A is an affordable wheelchair with 3 different styles to meet different customer's needs.


With a new shaped frame, which is completely different from traditional wheelchair. Height and Angel adjustable headrest offers better comfort for passengers.

Features

1. Reclining backrest with a simple button.
2. Detachable Headrest increases the comfort of the chair.
3. Brushless motor 180W, free wheel mode available.
4. Flip up armrest and swing away legrest.
5. With imported Lithium battery.
6. New designed wheel with absorbing ability.

Colors

Cushion color

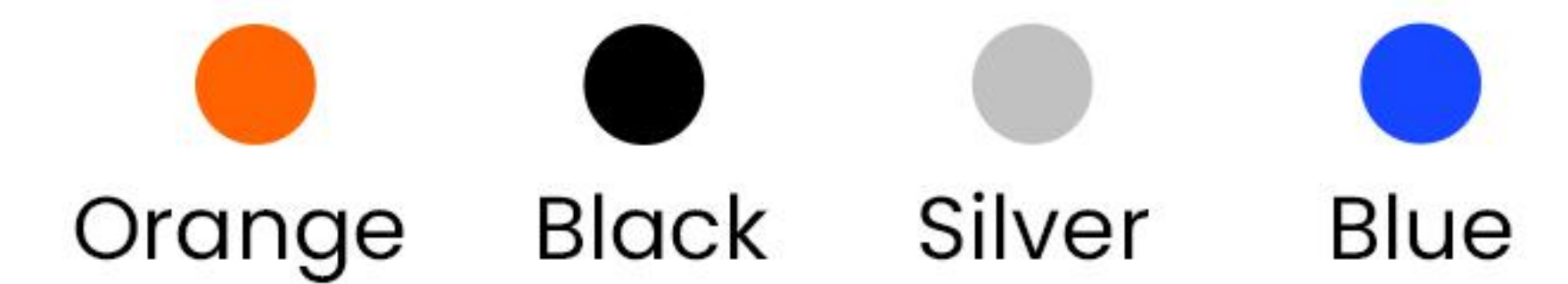

Frame color




Reclining steps



- Pull apart the left and the right sides buckles.
- Gravity makes passenger's back push down the backrest.
- Release the buckles till the backrest to an ideal angle



Detachable legrest

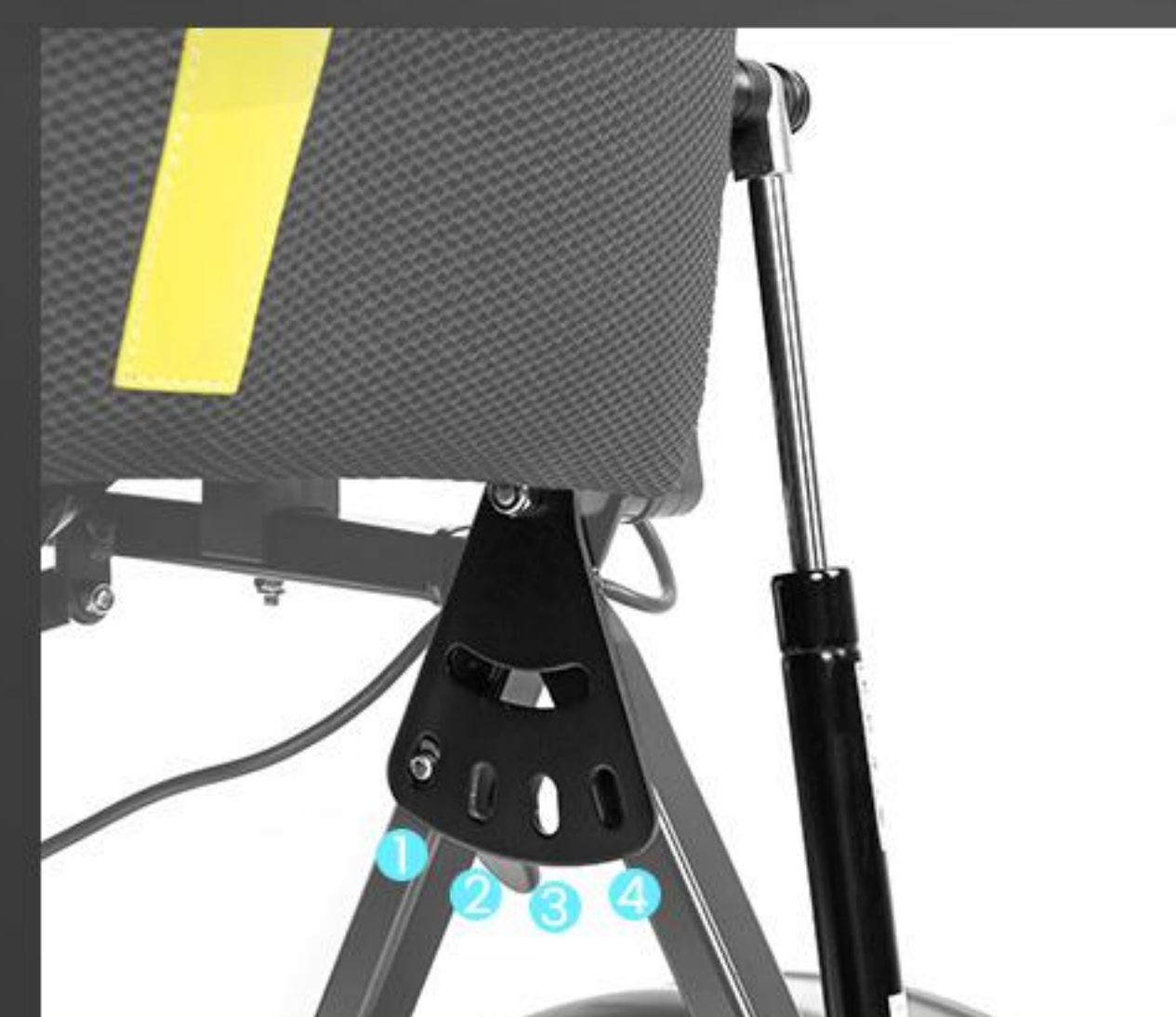

4 positions reclining adjustment



Module



Brushless motor



• YATTL INDUSTRY CO., LTD.

Tel: +86 755 8170 7199

Email: info@yattll.com